From Eye to Insight

0



OPTIMIZE YOUR WORKFLOW IN INSPECTION & REWORK

S9 i, S9 D, and S9 E stereo microscopes

000

DEVELOPED TO BOOST EFFICIENCY

S9

C

0

OPTIMIZE YOUR WORKFLOW

Work up to 20% faster with S9 series stereo microscopes*

Continuously improving production, keeping defect rates low, and fulfilling customer requests in order to stay competitive can be very challenging. Leica has developed the S9 stereo microscope series to help you cope with these challenges. The microscopes feature outstanding optical quality and FusionOptics technology for a three times greater depth of field. With a larger area in sharp focus, operators can immediately identify defects with fewer time-consuming microscope adjustments.

Eliminate extra steps in your workflow...

Locate area of interest

Take time to focus and zoom-in to see the details Go deeper, refocus and zoom-in again

Access the sample for reworks Document and transfer the results

Example of a standard workflow in stereo microscope inspections*

... boost efficiency with the S9 series and save up to 20% of your time*



EVERYTHING AT A GLANCE

See details faster to locate issues


No refocusing is needed to see all details within a height of up to 12 mm sharply in focus. The unique FusionOptics technology from Leica overcomes optical limitations. Operators will experience a more natural viewing perception, allowing them to work more efficiently.

Shift quickly from overview to detail

The S series offers high magnification of up to 55x in combination with a large 9:1 zoom range. This limits workflow interruptions from microscope adjustments and enables operators to inspect parts of different sizes more quickly.

FusionOptics technology from Leica

Conventional stereo microscopes have two identical beam paths that reveal a spatial impression of the sample. FusionOptics technology takes advantage of a neurological phenomenon: the microscope's **left beam path shows** an image with **great depth of field**, while the **right beam path** shows an image with **high resolution**. The human **brain** then effortlessly **combines the information from both channels into one image**. This results in an image perception with high resolution and a great depth of field at the same time – an unparalleled Leica standard in stereo microscopy.

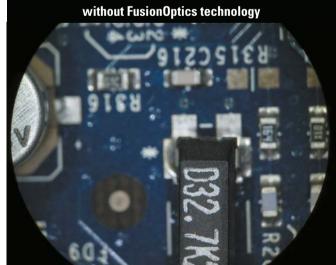

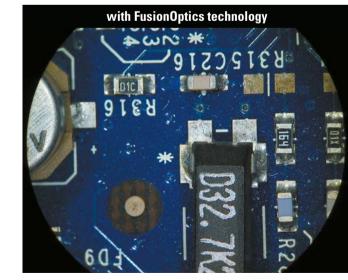


Image simulation, showing a printed circuit board sample without and with the FusionOptics effect, perceived when looking through the oculars.

3D VIEW WITH 12 MM DEPTH OF FIELD 37.7 MM OBJECT FIELD

SWIFTLY ACCESS YOUR SAMPLE

Maximum working distance for maximum speed

Sample inspection and manipulation under the microscope can be conveniently performed thanks to the 122 mm working distance of all S9 stereo microscopes. Enjoy working with more space for your tools under the microscope lens – sometimes every mm counts. Operators can access the sample with ease when using standardized tweezers or other manipulation tools.



Customize it and make it yours

In combination with different stands the S9 microscopes can be turned into customized workstations. This helps you to increase efficiency and the quality of your work results. The 35° viewing angle of the microscope allows a natural head posture and helps to prevent operator fatigue. For recurring tasks the switchable zoom click-stops assure reproducible and reliable results from operator to operator.

S9 stereo microscopes are ideally suited for production and inspection tasks for a variety of samples. The swingarm and flexarm offer even greater flexibility for different application areas, such as:

- > Electronic parts
- > Precision mechanics
- > Automotive engineering
- > Plastic parts
- > Forensic investigations
- > Medical device manufacturing
- Specimen dissection & preparation in life sciences









Flexarm stand for large samples

SAVE TIME: SHARE RESULTS DIGITALLY


Instant sharing for instant feedback

The S9 i stereo microscope with integrated 10 MP camera enables you to view images digitally at up to 35 frames per second on a PC, HD-monitor, or mobile device. This allows you to quickly and easily react to queries, get a second opinion, and discuss problems with colleagues. To work with mobile devices download the AirLab App from iTunes or the Google Play Store. For Apple Mac download Acquire.


Manage your images with ease


Inspect, analyze, and organize your images with the LAS X software from Leica. This latest software platform enables straightforward imaging and documentation, available for industrial and life science applications. It supports operators to deliver reliable results with confidence. LAS X also offers a range of additional software modules and expert solutions for specific applications.



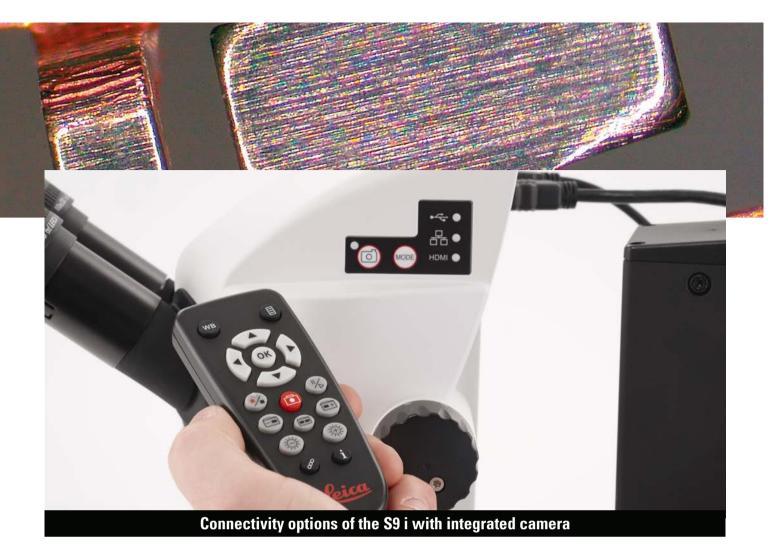


Use **USB mode** to connect the camera directly with a **PC** or Mac **via USB2** cable. Use **Ethernet mode** to connect the camera **with your facility's network** via Ethernet (RJ45) cable. Use **HDMI mode** to connect the camera to a large HD-monitor for a **standalone operation** without computer.

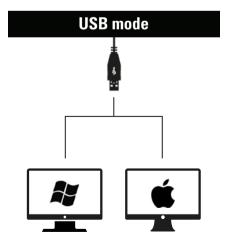

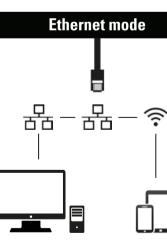

Capture images from any remote PC or Mac with access to your LAN network. Mobile devices can also be used, if they have access to your network ia Wi-Fi.

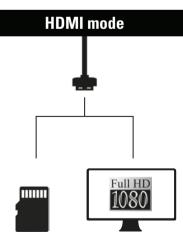

THE S9 SERIES

Viewing only

Inspect, observe, or manipulate your sample exclusively through the eyepieces. The **S9E** provides a cost effective solution with fast return on investment if you don't need to document with your microscope.


Sharing digitally


Share, document, and report results quickly and reliably. The integrated 10 MP CMOS camera of the **S9 i** can livestream images via USB, HDMI or Ethernet connection to your PC, HD-monitor, or mobile device.

Documentation ready

Stay flexible – add a camera to the **S9D** for documentation and sharing at any time. Observe your sample through the eyepieces and capture images simultaneously.



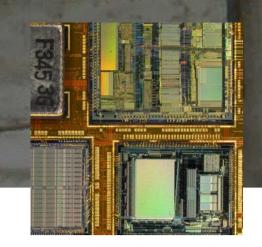

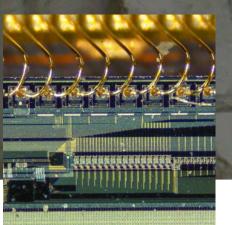

Prepare for the ongoing miniaturization trend in Electronics with 55x magnification and 9:1 zoom range to quickly change from overview to detail.


Put your Automotive sample in the right light with a wide range of illumination accessories to reveal defects.









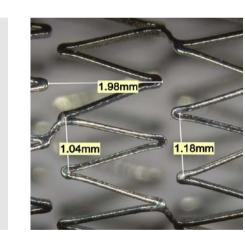





Sort and screen model organisms easily in Life Science Applications with the large working distance of 122 mm.



SPECIFICATIONS

Optical system, lead-free10° Greenough using best-corrected central part of the objective; complete apochromatically corrected microscope systemZoom9:1, apochromaticViewing angle35°ESD protectionAntistaticSpecific surface resistivity2×10 ¹¹ Ω/ square, discharge time <2 seconds from 1,000 V to 100 V
Viewing angle35°ESD protectionAntistaticSpecific surface resistivity2×10 ¹¹ Ω/ square, discharge time <2 seconds from 1,000 V to 100 V
ESD protectionAntistaticSpecific surface resistivity2×10 ¹¹ Ω/square, discharge time <2 seconds from 1,000 V to 100 VMagnification range (basic outfit)6.1 x - 55 xMaximum resolution500 lp/mmMaximum numerical aperture0.167Working distance (basic outfit)122 mmObject field diameter37.7 mmAdjustable zoom limitsClick-stops 10 x, 20 x, 30 x, 40 x, and 50 x
Specific surface resistivity2×10 ¹¹ Ω/square, discharge time <2 seconds from 1,000 V to 100 V
Magnification range (basic outfit) $6.1 x - 55 x$ Maximum resolution $500 lp/mm$ Maximum numerical aperture 0.167 Working distance (basic outfit) $122 mm$ Object field diameter $37.7 mm$ Adjustable zoom limitsClick-stops 10 x, 20 x, 30 x, 40 x, and 50 x
(basic outfit)Maximum resolution500 lp/mmMaximum numerical aperture0.167Working distance (basic outfit)122 mmObject field diameter37.7 mmAdjustable zoom limitsClick-stops 10x, 20x, 30x, 40x, and 50x
Maximum numerical aperture0.167Working distance (basic outfit)122 mmObject field diameter37.7 mmAdjustable zoom limitsClick-stops 10 x, 20 x, 30 x, 40 x, and 50 x
apertureWorking distance (basic outfit)0bject field diameter37.7 mmAdjustable zoom limitsClick-stops 10 x, 20 x, 30 x, 40 x, and 50 x
(basic outfit) Object field diameter 37.7 mm Adjustable zoom limits Click-stops 10x, 20x, 30x, 40x, and 50x
Adjustable zoom limits Click-stops 10 x, 20 x, 30 x, 40 x, and 50 x
Video/photo outlet,-50 % video-switchable50 % visual,permanent
Integrated camera-10 MP resolutionLive image up to 35 fps (1,024 x 768 pixels)Live image up to 35 fps (1,024 x 768 pixels)Sensor size 6.44 mm x 4.6 mm, 1/2.3" CMOSPixel size 1.67 µm x 1.67 µm
Standard objectives, lead-free Apochromats 0.5 x, 0.63 x, 0.75 x, 1.6 x, 2.0 x
Ergonomic eyepieces, 10x/23, 16x/16, 20x/12 fixed and adjustable, with cups
Ergonomic eyepieces10x/23, 16x/15,for eyeglass wearers,25x/9.5, 40x/6adjustable, with eyecups25x/9.5, 40x/6
Interpupillary distance 50 – 76 mm

Leica Microsystems (Schweiz) AG \cdot Max-Schmidheiny-Strasse 201 \cdot 9435 Heerbrugg, Switzerland T +41 71 726 34 34 \cdot F +41 71 726 34 44



CONNECT WITH US!

www.leica-microsystems.com